

**ASAP**<sup>®</sup>

**ASAPLP**<sup>™</sup>

Aspiration Catheter Kits

# POWERFUL CLOT ASPIRATION




## Powerful Clot Aspiration

ASAP® and ASAPLP™ Aspiration Catheter Kits are designed for the efficient removal of fresh, soft thrombi in PCI patients with large residual thrombus burden. Tailored to streamline procedures, ASAP and ASAPLP Aspiration Catheter Kits provide quality products you need to deliver powerful clot aspiration, restoring and improving blood flow in thrombosed vessels!

ASAP Aspiration Catheter	
Aspiration Lumen Area	1.01 mm <sup>2</sup>
Max Guide Wire Size	0.014", 0.36 mm
Catheter Length	140 cm
Catheter OD	0.068", 1.73 mm
Guide Catheter Compatibility	≥0.070", 1.78 mm, 6F

ASAPLP Aspiration Catheter	
Aspiration Lumen Area	0.81 mm <sup>2</sup>
Max Guide Wire Size	0.014", 0.36 mm
Catheter Length	145 cm
Catheter OD	0.055", 1.40 mm
Guide Catheter Compatibility	≥0.066", 1.68 mm, 6F

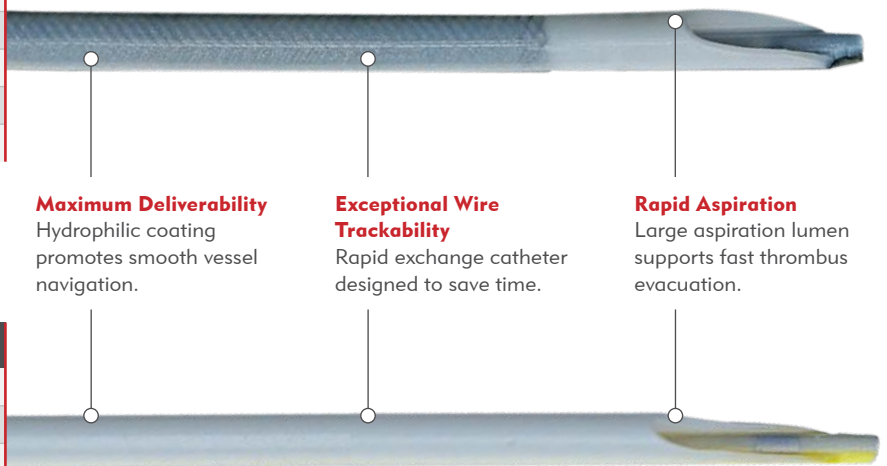

**Maximum Deliverability**  
Hydrophilic coating promotes smooth vessel navigation.



**Exceptional Wire Trackability**  
Rapid exchange catheter designed to save time.

**Rapid Aspiration**  
Large aspiration lumen supports fast thrombus evacuation.

### MORE KEY PRODUCT FEATURES & BENEFITS


- **Enhanced Angiographic Visualization** Depth markers and a radiopaque tip facilitate target vessel treatment.
- **Smooth Transitions** Tapered tip helps to maneuver through blockages.

Built with 100% stainless steel, the ASAP and ASAPLP Aspiration Catheters' braided technology breaks through clots without kinking, facilitating improved blood flow (TIMI ≥2).<sup>1</sup>

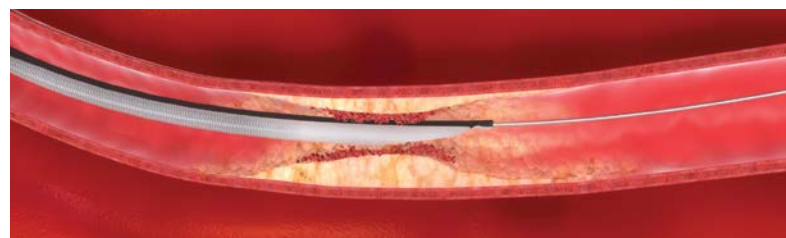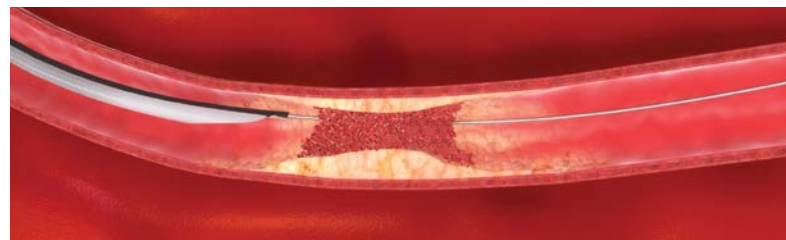


## PROVEN RESULTS IN A CONVENIENT KIT

Remove clots safely, quickly, and effectively with ASAP and ASAPLP Aspiration Catheters.

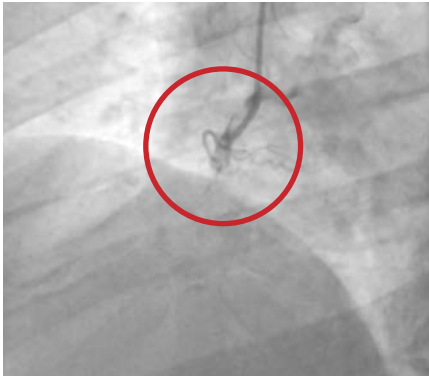

Large Thrombus Burden, Pre Aspiration

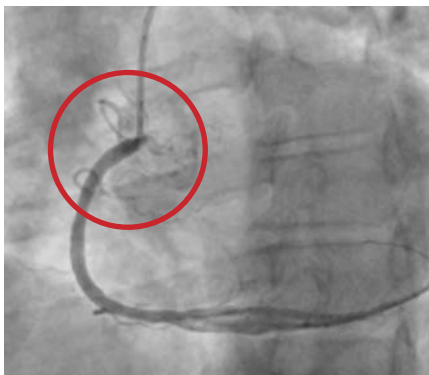

No Thrombus, Post Aspiration with ASAP



### ASAP Aspiration Catheter Kit

- 1- ASAP Aspiration Catheter
- 1- Extension Tubing Set with One-Way Stopcock, 21.5 cm total length
- 1- RXP® Flushing Syringe, 4 mL
- 1- MicroStop Fluid Collection Basin with Lid
- 2- Pore Filter Baskets, 70 micron
- 2- VacLok® Vacuum Pressure Syringes, 30 mL



### ASAPLP Aspiration Catheter Kit

- 1- ASAPLP Aspiration Catheter
- 1- Stiffening Stylet, preloaded in aspiration lumen of catheter
- 1- Extension Tubing Set with One-Way Stopcock, 21.5 cm total length
- 1- RXP® Flushing Syringe, 4 mL
- 1- MicroStop Fluid Collection Basin with Lid
- 2- Pore Filter Baskets, 70 micron
- 2- VacLok® Vacuum Pressure Syringes, 30 mL

## ORDERING INFORMATION

Catalog Number	Product Description	Units Per Box
ASAP100	ASAP Aspiration Catheter & Kit	1
ASAPLP	ASAPLP Aspiration Catheter & Kit	1

# COMBINE MERIT PRODUCTS. CREATE PROCEDURAL SOLUTIONS.


Choose from a wide selection of Merit products that are compatible with the ASAP and ASAPLP Aspiration Catheter Kits. Each works to provide you with the complete procedural solutions you need to achieve better patient outcomes.

## GUIDE WIRES AND ACCESSORIES


**Merit Laureate®**  
Hydrophilic  
Guide Wires



**Tenor®**  
Steerable  
Guidewires  
0.014" (165 cm)

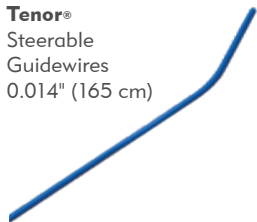


**InQwire®**  
Diagnostic Guide Wires

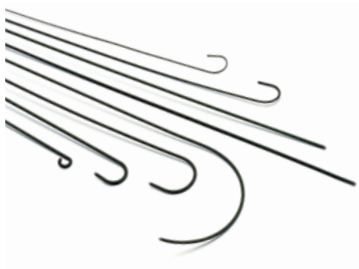

## GUIDING CATHETERS

**ConcierGE®**  
Guiding Catheters

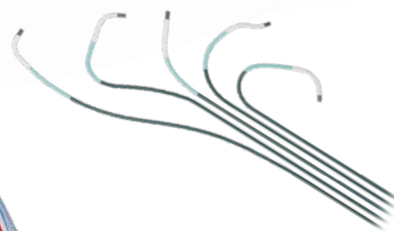

## ANGIOGRAPHY CATHETERS

**Performa®**  
Diagnostic Cardiology  
Catheter



## ACCESS

**Prelude IDEAL™**  
Hydrophilic Sheath  
Introducer



1. Data on File.

Before using, refer to Instructions for Use for indications, contraindications, warnings, precautions, and directions for use.



Understand. Innovate. Deliver.™

[merit.com](http://merit.com)

Merit Medical Systems, Inc.  
1600 West Merit Parkway  
South Jordan, Utah 84095  
Main: +1.801.253.1600  
Customer Service: +1.801.208.4300

Merit Medical Europe, Middle  
East, & Africa (EMEA)  
Amerikalaan 42, 6199 AE  
Maastricht-Airport  
The Netherlands  
+31 43 358 82 22

Merit Medical Ireland Ltd.  
Parkmore Business Park West  
Galway, Ireland  
+353 (0) 91 703 733